


# MIRUS™ Ligation | Titanium

The MIRUS™ Ligation Clips are implantable medical devices. They are widely used in surgeries for hemostasis as well as for marking purposes. They are designed for use on tubular structures or blood vessels, wherever indicated. Ligation clips are commonly used in cardiothoracic, vascular, gastrointestinal, urological, gynecological, bariatric and general procedures.

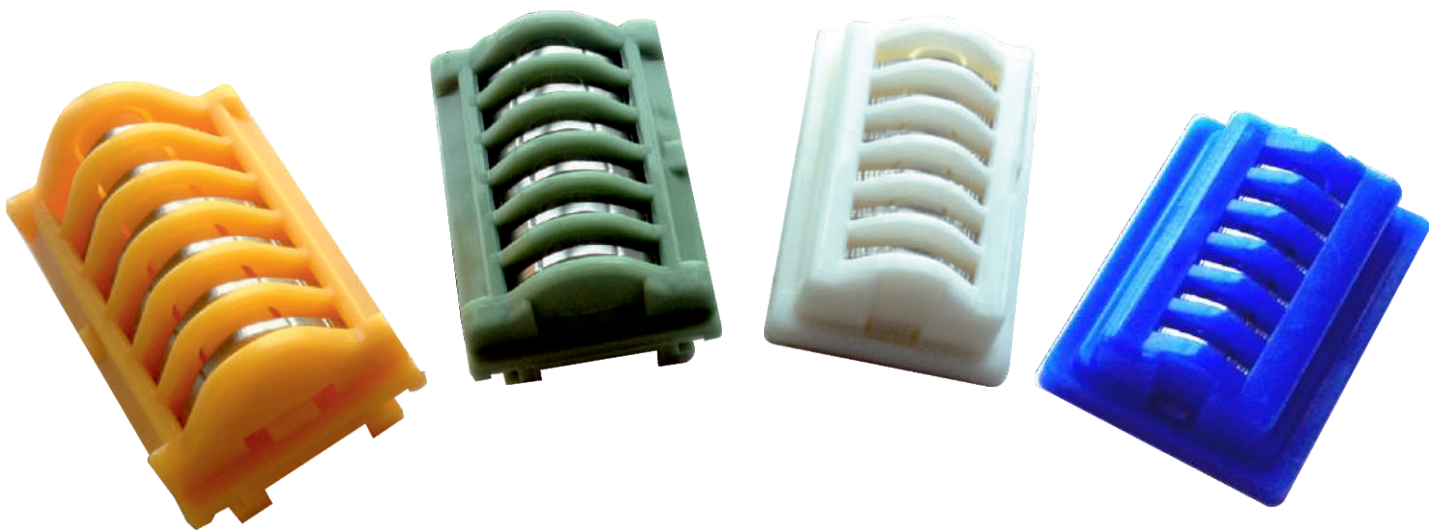

## Meril Titanium Clips:

The titanium clip engulfs tissue with precise tip-to-tip closure. The clip wire is designed to give each clip a firm grip on vessels, while the triangulated cross-section of the clip leg maximizes surface-to-surface contact between clip and jaw, practically eliminating clip fallout.

Lateral and Transverse grooves allow secure fixation on the structure and increase resistance to dislodgement of a formed clip.

MIRUS™ ligation clips can be used in endoscopic procedures as well as in open surgeries. They are compatible with Meril applicator and other applicators available in the market.

## Available as :

Product code	Clip size	Clips/Cartridge	Cartridges/Box	Clips/Box
MLT100	Small	6	20	120
MLT200	Medium	6	20	120
MLT300	Medium-Large	6	20	120
MLT400	Large	6	20	120

Ligation Clips

Titanium

Procedure Application

Blood vessel ligation for hemostatic purposes

# MIRUS™ Ligation | Titanium

## Appliers for Open Surgery:



Product	Product code	Specification	Length	Applicator size	Clips compatibility
Applicator-Open	FTOA-10015	OPEN APPLICATOR 100-15CM	15 cm	Small	MLT100
Applicator-Open	FTOA-10020	OPEN APPLICATOR 100-20CM	20 cm	Small	MLT100
Applicator-Open	FTOA-10028	OPEN APPLICATOR 100-28CM	28 cm	Small	MLT100
Applicator-Open	FTOA-20015	OPEN APPLICATOR 200-15CM	15 cm	Medium	MLT200
Applicator-Open	FTOA-20020	OPEN APPLICATOR 200-20CM	20 cm	Medium	MLT200
Applicator-Open	FTOA-20028	OPEN APPLICATOR 200-28CM	28 cm	Medium	MLT200
Applicator-Open	FTOA-30015	OPEN APPLICATOR 300-15CM	15 cm	Medium-Large	MLT300
Applicator-Open	FTOA-30020	OPEN APPLICATOR 300-20CM	20 cm	Medium-Large	MLT300
Applicator-Open	FTOA-30028	OPEN APPLICATOR 300-28CM	28 cm	Medium-Large	MLT300
Applicator-Open	FTOA-40015	OPEN APPLICATOR 400-15CM	15 cm	Large	MLT400
Applicator-Open	FTOA-40020	OPEN APPLICATOR 400-20CM	20 cm	Large	MLT400
Applicator-Open	FTOA-40028	OPEN APPLICATOR 400-28CM	28 cm	Large	MLT400

## Appliers for Laparoscopic Surgery:



Product	Product code	Specification	Length	Applicator size	Clips compatibility
Applicator-Endo	FTEA-00300	LAPAROSCOPIC APPLICATOR 300	Standard	Medium-Large	MLT300
Applicator-Endo	FTEA-00400	LAPAROSCOPIC APPLICATOR 400	Standard	Large	MLT400

# MIRUS™ Ligation | Titanium

## V-Shape Ligation Clip

VMLT-060



VMLT-080



VMLT-100




VMLT-200




VMLT-300



VMLT-400

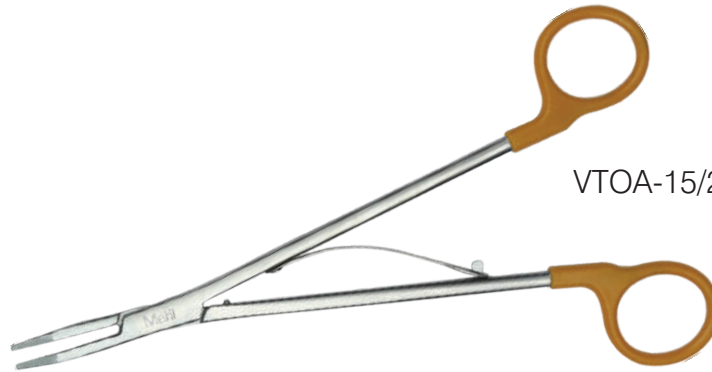

VTOA-15/20/28

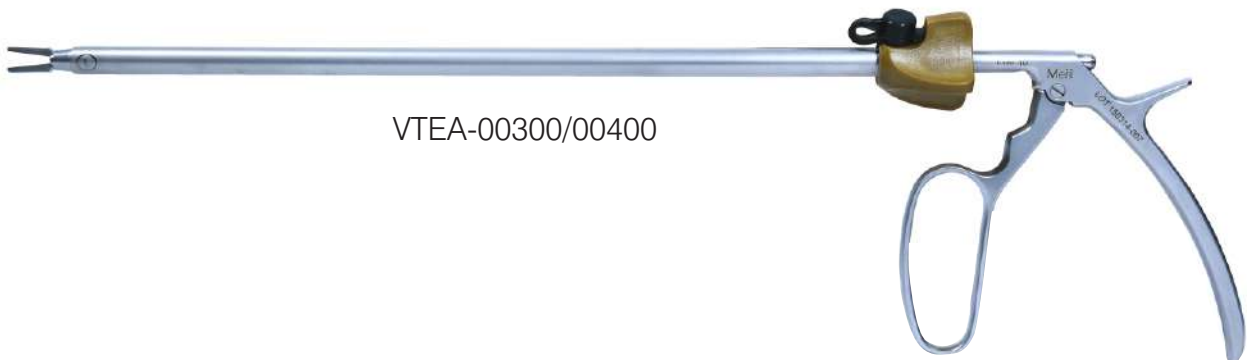

VTEA-00300/00400

Product Code	Clip Size	Clip/Cartridges	Cartridges /Box	Clips/Box
VMLT-060	Micro	6	20	120
VMLT-080	Small	6	20	120
VMLT-100	Small-Medium	6	20	120
VMLT-200	Medium	6	20	120
VMLT-300	Medium-Large	6	20	120
VMLT-400	Large	6	20	120

Ligation Clips	Titanium
Procedure Application	Blood vessel ligation for hemostatic purposes